



# sonicsmooth

by MICHAEL TODD BEAUTY

Sonic Dermaplaning,  
Exfoliation & Facial Hair  
Removal System

INSTRUCTIONS  
FOR USE

Digital Owner's  
Manual





# CONGRATULATIONS

Dear Customer,

Thank you and congratulations on choosing Michael Todd Beauty's SonicSmooth, one of the safest and most effective dermaplaning systems available for home use. We are confident that you will adore your instantly smoother, brighter and "fuzz free" skin for years to come.

SonicSmooth is the world's first sonic dermaplaning system specifically designed to expertly trim facial vellus hair (peach fuzz) while helping to gently remove dead skin cells and smooth away rough, uneven textured skin.

The SonicSmooth sonic dermaplaning system uses advanced sonic wave technology paired with Michael Todd Beauty's dermaplaning blade to expertly provide safe, instant and painless professional-grade dermaplaning results in the comfort of your own home.

3 Speed settings allow you dial up or dial down the power/intensity setting on the SonicSmooth so that you are able to customize your dermaplaning session to your specific speed preferences.

SonicSmooth is cordless and rechargeable with a 3 hour max. running battery life. The device handle is coated in a soft rubber for comfort, to prevent slipping and for ease of use. SonicSmooth has an ergonomically designed handle that comfortably sits securely in the palm of your hand.

SonicSmooth will quickly become your favorite weekly skin smoothing treatment. You are going to love the silky soft feeling of "fuzz free" skin and will immediately notice how much brighter your complexion looks. Your serums and creams will be able to work better after a SonicSmooth session because the dead skin cells and build-up will be gently sloughed away by state-of-the-art sonic technology. In addition, your makeup is going to sit flawlessly on your skin after removing the peach fuzz and layers of dead skin. ENJOY!

READ ALL INSTRUCTIONS BEFORE USE  
AND SAVE FOR FUTURE REFERENCE

## REGULATORY INFORMATION

U.S. This Class B digital product complies with Part 15 of the U.S. Federal Communications Commission (FCC) rules. Operation is subject to the following two conditions: (1) this device may not cause interference, and (2) this device must accept any interference, including interference that may cause undesired operation of the device.

Canada. This device complies with Industry Canada licence-exempt RSS standard(s). Operation is subject to the following two conditions: (1) this device may not cause interference and (2) this device must accept any interference including interference that may cause undesired operation of the device.

Le présent appareil est conforme aux CNR d'Industrie Canada applicables aux appareils radio exempts de licence. L'exploitation est autorisée aux deux conditions suivantes : (1) l'appareil ne doit pas produire de brouillage, et (2) l'utilisateur de l'appareil doit accepter tout brouillage radioélectrique subi, même si le brouillage est susceptible d'en compromettre le fonctionnement."

Mfg: Michael Todd Beauty

Model:  
Sonicsmooth  
Device complies with  
CAN ICES-3 (B)/NMB-3(B)

For home or office use. Not intended for use in machinery, medical or industrial applications.

Radio and TV Interference Regulations.  
This product has been tested and found to comply with the limits for a Class B digital device pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This product generates, uses and can radiate radio frequency (RF) energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications devices (such as televisions, baby monitors, cordless phones, etc). However, there is no guarantee that RF interference will not occur in a particular installation. If this product does cause harmful interference to radio or television reception, which can be determined by turning the product off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antennae of the other radio-communications device
- Move the equipment farther away from the radio or tv
- Connect the product into an outlet on a circuit different from that to which the other device is connected
- Consult the manufacturer's customer service or an experienced radio/tv technician for more suggestions.

# IMPORTANT SAFETY INFORMATION

BASIC SAFETY PRECAUTIONS. THIS PRODUCT IS NOT DESIGNED OR INTENDED PRIMARILY FOR CHILDREN 12 YEARS OF AGE OR YOUNGER. When using electrical products, especially when children are present, basic safety precautions should always be taken. IF YOU THINK THE POWER IS OFF WHEN THE SWITCH IS OFF, YOU ARE WRONG. If the plug is in, the power is on. Unplug electrical appliances when you are not using them.

## DANGERS. To reduce the risk of electrocution:

1. **DO NOT** place or store the Charging Cradle where it can fall or be pulled into a tub or sink.
2. **DO NOT** place the charging cradle or plug into water or other liquid. The charging cradle should not be immersed or used in shower or bath.
3. **DO NOT** reach for the Charging Cradle or the plug that has fallen into water. If unplug while in water, might get electrocuted.
4. **DO NOT** use Sonicsmooth in a location where the power cord can come in contact with water while it is plugged in.
5. Keep Sonicsmooth away from open flame or flammable aerosol products. Sonicsmooth is not suitable for use in the presence of a flammable anesthetic mixture with air or oxygen or nitrous oxide.

**WARNINGS. To reduce the risk of burns, electrocution, fire or injury to persons:**

1. Close supervision is necessary when this product is used by, on, or near children or physically challenged persons.
2. Use this product only for its intended use as described in these Instructions for Use. **DO NOT** use attachments other than those recommended and provided by MICHAEL TODD BEAUTY.
3. Keep the power cord away from heated surfaces.
4. Never use Sonicsmooth while sleeping or drowsy.
5. **DO NOT** use Sonicsmooth outdoors or where oxygen is being administered.
6. Never operate Sonicsmooth if it is damaged.
7. **DO NOT** use the charging cradle for recharging if the power cord or plug is damaged, the charging cradle is not working properly, or if the charging cradle has been dropped on the ground or into water.
8. The handle and charging cradle have no serviceable parts. Never attempt to open the handle for any reason. If the charging cradle cord is damaged, it must be unplugged and replaced. Contact customer service at +1-772.446.0149 or email us at [customerservice@michaeltoddbeauty.com](mailto:customerservice@michaeltoddbeauty.com) for instructions on how to obtain a replacement.

9. Use Sonicsmooth with the supplied power adapter only.

#### IMPORTANT CAUTIONS.

1. If you experience any pain or discomfort from using Sonicsmooth, discontinue use.
2. To avoid risk of injury, we **DO NOT** recommend using Sonicsmooth over the eyelid, directly on the nose or on the lips.

# CONTENTS

Your Sonicsmooth comes with the following:

- A. Sonicsmooth Dermaplaning Tool
- B. Stainless Steel, Single-Use Dermaplaning Blades - (8wks supply)
- C. USB Cable
- D. Charging-Storage Cradle



# GETTING TO KNOW YOUR SONICSMOOTH



# OPERATING INSTRUCTIONS

## 1. Prep the Skin

Thoroughly cleanse the treatment area with an oil-free cleanser, then rinse and pat dry. Sonicsmooth is most effective on clean, dry skin. If you have dry or sensitive skin, you can use a lightweight, non-comedogenic facial oil during dermaplaning. This creates a barrier that reduces friction and irritation while enhancing hydration, making the process gentler and more comfortable.

## 2. Charge the Device

Before using your Sonicsmooth for the first time, charge the handle for 12 hours. Refer to the charging instructions for details.

## 3. Areas to Avoid

Do not use the Sonicsmooth directly on the eyelids, nose, or lips (upper lip and below the lips are fine for treatment).

## 4. Attaching and Removing Dermaplaning Blades:

- Always ensure the device is OFF when not in use.
- Use the cartridge case to attach and remove dermaplaning blades:
  - **To attach a new dermaplaning blade:**  
Slide the device handle along the tracks of a new blade until it clicks into place.
  - **To remove a used blade:**  
Insert the blade into the cartridge case until you hear a click, then slide the device handle away from the case.
- Never use your hands or fingers to handle blades.

- ## 5. Once the blade is securely attached, lift the device with the attached blade from the cartridge case, and you're ready to begin your treatment.



6. Starting: Press the ON/OFF Button to turn on
7. Selecting Speeds: To customize your dermaplaning session, Sonicsmooth provides a choice of 3 power/intensity settings (low, medium and high). We recommend that you start on a lower speed setting and gradually increase power/intensity as you get used to sonic dermaplaning.


After pressing the on/off button one time the unit will power on in low speed. When you press the button again, the unit will change to medium speed. Pressing the power button a third time will change the speed to high and pressing the button one more time will shut off the device.

Caution: Pull hair back away from face and take precaution when using around hairline and eyebrows.

Sonicsmooth is for external use only.

Do not use on broken or irritated skin.

This device is for cosmetic use only and not intended for medical use.

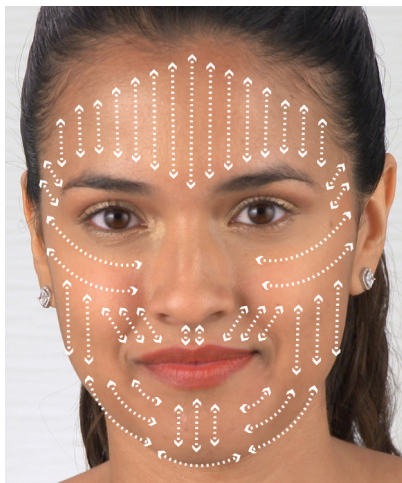

## Sonicsmooth Dermaplaning Tips

This is a graphic that shows the aesthetician recommended directions to dermaplane in different areas of the face for the best possible result.

8. Starting at the top of your face near the forehead hairline, begin making short feathered strokes following the direction of the hair growth. You may gently glide over skin in a treatment area 1-3 times using the blade at a 45° angle. We DO NOT recommend using Sonicsmooth over raised imperfections (including acne) the eyelid, directly on the nose or on the lips.
9. When you are finished, power off the Sonicsmooth, discard your dermaplaning blade and return it to the charging cradle.
10. After dermaplaning your skin may be slightly pink or red in color. This is normal. After treatment, you may use serums or moisturizers. DO NOT use scrubs, exfoliators or peels. Any sensation of post treatment redness will fade in the next couple of hours.

## Expert Tips from an Aesthetician

1. It is important to space treatments apart every 7-10 days.
2. It's important to protect your skin from sun exposure, harsh chemicals and exfoliants around your dermaplaning sessions.
3. Ideally, performing SonicSmooth treatments at night before bed will allow your skin to refresh itself while you sleep. Keep in mind that after dermaplaning skin is more readily able to absorb serums and moisturizers because the dead layer of skin and vellus hair that blocks absorption has been safely removed. You may notice a dramatic improvement in the look and feel of your skin when your serums and creams are able to perform properly.
4. It is not recommended that you apply makeup directly after a SonicSmooth session. Liquid, cream and powder makeup will have a smooth fresh palette elevating the appearance of makeup when applied the morning after or hours after a SonicSmooth session.
5. We do not recommend using the SonicSmooth in the eye socket area, on the nose or directly on the lips. These areas are not flat in surface and are not suitable for the flat dermaplaning blade on the SonicSmooth.

6. Holding your Sonicsmooth at a 45° angle is the safest and most effective way to sonic dermaplane. It is not recommended that you hold the Sonicsmooth perpendicular to skin.
7. Unlike shaving other parts of your body, the best possible result from Sonicsmooth is achieved by using light feathery strokes in a downward motion which is the same direction as hair growth. When we shave other parts of our bodies we typically shave against the direction of hair growth but when dealing with soft vellus facial hair, dermaplaning in the same direction as hair growth yields a better result.
8. When working in areas better accessed in a slightly sideways angle, it is suitable to angle the Sonicsmooth to fit the contour of the face. An example of this is when working on and around the cheek bone area. It may be easier to successfully dermaplane this area by directing strokes from the hairline towards the center of the face. Just remember to hold skin taut with the opposite hand.


# CHARGING

1. Your Sonicsmooth is powered by a rechargeable battery. Sonicsmooth should be recharged after each use. Depending on the power setting used, the fully charged Sonicsmooth should provide 3 hours of continuous use before needing to be recharged.
2. First Use. Before you use your Sonicsmooth for the first time, you need to charge it fully for 12 hours. For subsequent recharging allow 8-10 hours.
3. Attach the included USB Charging Cable into a USB port and connect to a 110v wall adapter. Place Sonicsmooth in the Charging Cradle with the Power Button facing forward. The Michael Todd Beauty logo on the Charging Cradle and the Power Button should be facing the same direction.
4. The Charge Indicator will light up in red indicating the unit is charging. When the unit is finished charging the indicator will light up in green. Unplug the adapter from the wall and Charging Cradle when not in use. DO NOT unplug the adapter from the wall outlet when it is still connected to the unit as this will cause discharging of the battery.

5. Infrequent Use. The Handle needs to receive a charge occasionally. If you rarely use the Handle, recharge it in the Charging Cradle every 3 months to maintain its charging capability.
6. Charge Indicator. When the charge is good the Power Light will be green. When the charge is low the Power Light will be red. After 10 minutes of use with a low charge the Power Light will Flash Red and the device will power off.

If the Sonicsmooth is not sitting in the Charging Cradle Properly it will fall out and will not receive a charge.

If your Sonicsmooth stops working, charge it in the Charging Cradle for a minimum of 5 hours. If the issue is not resolved, please contact Customer Service.

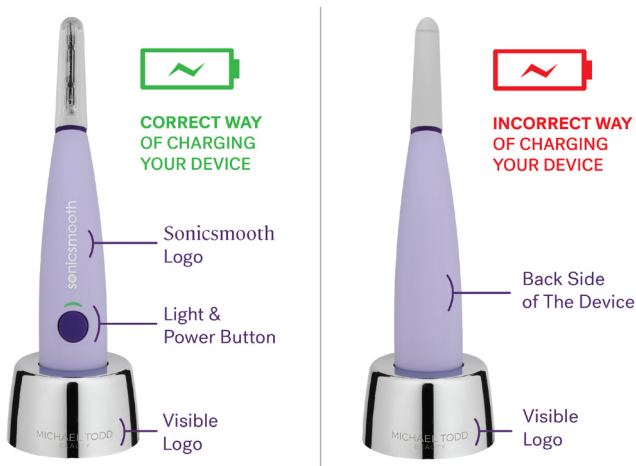

# CARE AND MAINTENANCE

1. Handle. To periodically clean the Handle, use an alcohol wipe. Using water on the handle may damage the unit.
2. The dermaplaning blades are designed for single use. Dermaplaning with a fresh dermaplaning blade every time is a good hygienic practice as well as eliminates issues associated with using a dull blade.
3. Charging Cradle. **DO NOT** immerse the Charging Cradle in water. Unplug the Charging Cradle before cleaning it. After unplugging the Charging Cradle, wipe with a soft, damp cloth.
4. **DO NOT** clean Sonicsmooth or any of its components using a dishwashing machine.
5. **DO NOT** use abrasive cleaners or chemicals to clean your Sonicsmooth.
6. Storage. When not in use for prolonged periods, store your Sonicsmooth in a dry place at a temperature between 32°F and 140°F (0°C to 60°C).
7. Sonicsmooth does not require any other regular maintenance. If the unit is not working properly, discontinue use and contact Customer Service at [customerservice@michaeltoddbeauty.com](mailto:customerservice@michaeltoddbeauty.com) or call +1-772.446.0149.

## Disposal Instructions

The rechargeable battery and device must not be disposed of in household waste. Dispose of in accordance with local regulations for electronic waste and battery recycling. Contact your local council or authorized recycling facility for guidance.

# LIMITED WARRANTY

Please use this as your reference for any customer service or warranty needs you may have in the future.

You must register your product within 90 days of purchase in order to obtain the warranty. Your warranty is not effective unless you register your product within the time provided. In order to register your product go on-line to [www.warrantymichaeltoddbeauty.com](http://www.warrantymichaeltoddbeauty.com) and complete the information required for registration.

LIMITED 3 YEAR WARRANTY. This limited warranty gives you specific legal rights and you may have other rights which vary from jurisdiction to jurisdiction. Any claim for breach of this limited warranty must be commenced within 3 years following the original purchase.

For a period of three years from the date of purchase, Michael Todd Beauty LP (“MICHAEL TODD BEAUTY”) warrants the Sonicsmooth (the “Product”) against Covered Defects. If there is a Covered Defect during the warranty period MICHAEL TODD BEAUTY agrees to exchange the Product with a new equivalent Product. This is the sole remedy available for breach of this limited warranty.

The limited warranty begins on the date of purchase and expires 3 years thereafter. In no event is the warranty period tolled or extended - including any period during which the Product (or any parts thereof) is in repair or if the Product is replaced by MICHAEL TODD BEAUTY during the warranty period.

This warranty is limited and applies only to a purchase of the Product from MICHAEL TODD BEAUTY or an authorized reseller. The warranty is for the original end user of the Product only and cannot be transferred under any circumstances.

Except as otherwise provided by law, to obtain service under this warranty, return the Product with a Covered Defect to us with your purchase receipt and, if you are returning the Product after one (1) year from the date of purchase, \$30 for postage and handling. For shipping instructions, contact our Customer Service department at [customerservice@michaeltoddbeauty.com](mailto:customerservice@michaeltoddbeauty.com) or call +1-772.446.0149.

WHAT IS COVERED? “Covered Defects” means damage, defects or failure resulting from ordinary consumer use under normal home conditions according to the instruction manual included with the Product. Covered Defects does NOT include using the Product with a voltage converter; damage, defects or failure caused by misuse or improper care of the Product; impact or insertion of foreign objects into or on the Product; improper storage; improper or incorrectly performed maintenance or repair of the Product; application of harmful chemicals or abrasives to the Product; failure to follow manufacturer’s instructions; normal wear and tear; alteration or overuse of the Product; as well as any professional or commercial use of the Product.

WARRANTY LIMITATIONS. THIS LIMITED WARRANTY IS EXCLUSIVE OF ALL OTHER WARRANTIES, WHETHER ORAL OR WRITTEN, EXPRESS OR IMPLIED. THERE ARE NO IMPLIED WARRANTIES CREATED BY THE MANUFACTURE, SALE OR USE OF THE PRODUCT.

IN NO EVENT SHALL MICHAEL TODD BEAUTY BE LIABLE FOR ANY SPECIAL, INCIDENTAL, CONSEQUENTIAL, PUNITIVE, OR INDIRECT DAMAGES BASED UPON BREACH OF WARRANTY, BREACH OF CONTRACT, NEGLIGENCE, STRICT PRODUCT LIABILITY, TORT OR ANY OTHER LEGAL THEORY.

THESE WARRANTY TERMS, CONDITIONS AND LIMITATIONS ALLOCATE THE RISKS OF PRODUCT FAILURE BETWEEN MICHAEL TODD BEAUTY AND YOU, THE END USER. THIS ALLOCATION IS RECOGNIZED BY BOTH PARTIES AND IS REFLECTED IN THE PRICE OF THE PRODUCT. BY USING SONICSMOOTH YOU ACKNOWLEDGE THAT YOU HAVE READ THESE TERMS, CONDITIONS AND LIMITATIONS, UNDERSTAND THEM AND ARE BOUND BY THEM.

SOME STATES DO NOT ALLOW LIMITS ON WARRANTIES, OR ON REMEDIES FOR BREACH IN CERTAIN TRANSACTIONS. IN SUCH STATES, THE LIMITS CONTAINED HEREIN AND THE DISCLAIMER OR EXCEPTIONS FROM WARRANTY MAY NOT APPLY. IF APPLICABLE LAW DOES NOT PERMIT THE DISCLAIMER OF IMPLIED WARRANTIES, THEN THE DURATION OF SUCH IMPLIED WARRANTIES, ARE EXPRESSLY LIMITED TO THE DURATION OF THE EXPRESS LIMITED WARRANTY SET FORTH HEREIN.

PRODUCT CHANGES. MICHAEL TODD BEAUTY reserves the right to make changes and improvements to its products without incurring any obligation to incorporate such improvements in products previously sold or to notify you of these changes or improvements.

DISPUTES. This limited warranty is governed by the laws of the State of Florida, without regard to principles of conflicts of laws, and applicable United States federal law. The United Nations Convention on Contracts for the International Sale of Goods shall not apply. YOU AGREE TO SUBMIT ANY CLAIM, DISPUTE, OR CONTROVERSY (REGARDLESS OF THE FORM OR TYPE OF ACTION) AGAINST US TO ARBITRATION ADMINISTERED BY THE NATIONAL ARBITRATION FORUM UNDER ITS CODE OF PROCEDURE THEN IN EFFECT (available at <http://www.arb-forum.com> or via telephone at 1-800-474-2371). Except to the extent otherwise required by law, you agree that the arbitration proceedings will take place exclusively where the principal office of MICHAEL TODD BEAUTY is located at the time of institution of the proceedings. You may be entitled to file suit under state laws without submitting your claim to arbitration first. Information on arbitration may be obtained from, and a claim may be filed with, the National Arbitration Forum at P.O. Box 50191, Minneapolis, MN 55405-0191 USA.



TRADE NAME: MICHAEL TODD BEAUTY

MODEL: SONICSMOOTH

Tested to Comply with FCC Standards



Michael Todd Beauty  
584 NW University Blvd, Suite 600  
Port St. Lucie, FL 34986, USA

Customerservice@Michaeltoddbeauty.com














Botley Mills, Southampton,  
SO30 2GB, UK

For more information, updates, reviews frequently asked questions and more visit us on the web at [www.michaeltoddbeauty.com](http://www.michaeltoddbeauty.com)

© 2026 MICHAEL TODD BEAUTY LP. All rights reserved.

# SYMBOLS EXPLAINED

Symbol	Meaning
	<b>Manufacturer</b> – Indicates the legal manufacturer of the device.
	<b>Caution / Warning</b> – Indicates that important safety information must be reviewed before use.
	<b>Consult Instructions for Use (IFU)</b> – Indicates the user should read the accompanying documentation before operating the device.
	<b>Regulatory Compliance Mark</b> – Indicates the product conforms to applicable regulatory requirements.
	<b>Keep Dry</b> – It indicates that the product should not be exposed to water or moisture during storage or transport.
	<b>CE Marking</b> – Indicates conformity with applicable European Union health, safety, and environmental protection legislation.
	<b>Authorized Representative in the United Kingdom</b> – Identifies the UK-based representative authorized to act on behalf of the manufacturer.
	<b>Batch Code / Lot Number</b> – Identifies the product's manufacturing batch for traceability purposes.
	<b>Federal Communications Commission Mark</b> – Indicates compliance with U.S. FCC requirements for electronic devices (radio frequency emissions).
	<b>Restriction of Hazardous Substances Compliance</b> – Indicates compliance with EU RoHS Directive limiting hazardous materials in electrical and electronic equipment.
	<b>WEEE Symbol (Waste Electrical and Electronic Equipment)</b> – Indicates the product must not be disposed of in household waste and must be collected separately for recycling.

